

Power, Versatility and Efficacy; One Beautiful Combination

To keep up with the ever-growing demand for aesthetic treatments, you need a system that offers a wide range of applications delivered effectively and efficiently. GentleMax gives you the power to offer today's most popular aesthetic treatments with the kind of superior performance you expect from an industry leader like Candela.

Superior Performance Means Greater Satisfaction

Efficacy counts. But so does efficiency. With GentleMax you can offer both. With the power of Candela technology, GentleMax can produce superior results in less time, and often with fewer treatments. Something that's sure to please patients and practitioners alike.

GentleFacial™

GentleFacial—a GentleMax Exclusive That Delivers Today's Most Popular Cosmetic Treatments

Wrinkles, facial veins, sun-damaged skin, uneven pigmentation, tone and texture, skin tightening—every patient has unique needs and expectations. With GentleMax, you can offer them a series of custom GentleFacial treatments designed around their individual needs.



See the Difference Only Candela Can Deliver.


One system—unlimited potential. With GentleMax, Candela has succeeded where others fail. It's simply what we do. For 37 years, we've led the way in laser and light-based technology—delivering not only superior products, but unmatched after-the-sale service and support. Customer service around the world is handled by our dedicated team of technicians. Our clinical support is recognized throughout the industry, and our Pathways purchasing program offers the choice, flexibility and value every practice will appreciate. And once you've made your purchase, MyCandela.com, our online marketing resource, gives you the tools to make the most of your Candela investment.

To learn more, visit www.candelalaser.com
or call 800-821-2013.

Optional Dynamic Cooling Device

Integrated controls, cryogen container and handpiece with distance gauge

Cryogen	HFC 134a
DCD Spray	User adjustable range:
Duration	10 – 100 ms
DCD Delay	User adjustable range:
Duration	3, 5, 10 – 150 ms
DCD Post	User adjustable range:
Spray Duration	0 – 50 ms

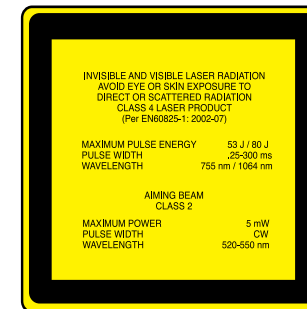

Complies with FDA performance standards for laser products except for deviations pursuant to Laser Notice No. 50 dated July 26, 2002.




Specifications:

Laser Type	Nd:YAG	Alexandrite
Wavelength	1064 nm	755 nm
Repetition	Up to 10 Hz	Up to 10 Hz
Pulse Duration	0.250 – 300 ms	0.250 – 300 ms
Spot Sizes (diameter)	1.5 mm, 3 mm, 6 mm, 8 mm, 10 mm, 12 mm, 15 mm, 18 mm	1.5 mm, 3 mm, 6 mm, 8 mm, 10 mm, 12 mm, 15 mm, 18 mm
Fluence Range	Up to 600 J/cm²	Up to 600 J/cm²
Beam Delivery	Lens-coupled optical fiber with handpiece	Lens-coupled optical fiber with handpiece
Pulse Control	Finger switch, Foot switch	Finger switch, Foot switch
Dimensions	45" H x 26.5" W x 38" D (114 cm x 67 cm x 97 cm)	
Weight	342 lbs (155.45 kgs)	
Electrical	230 V 50/60 Hz, single phase, 17.4A	

GentleYAG — Spot sizes 1.5 – 18 mm, repetition rates up to 10 Hz, 26,333 watts of peak power, fluences up to 600 J/cm². Microseconds to milliseconds pulse durations.

GentleLASE — Spot sizes 1.5 – 18 mm, repetition rates up to 10 Hz, 17,666 watts of peak power, fluence range up to 600 J/cm². Microseconds to milliseconds pulse durations.

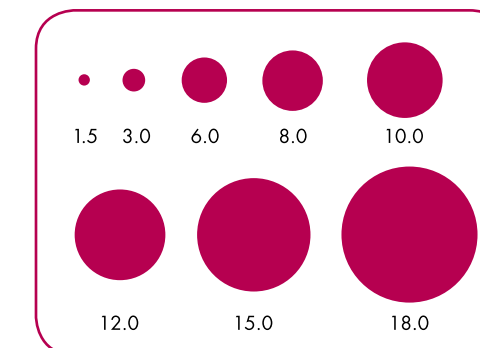

Candela Corporation

530 Boston Post Road
Wayland, MA 01778 USA
Phone: (508) 358-7637
Fax: (508) 358-5569
Toll Free: (800) 821-2013
www.candelalaser.com



GentleMax™



Integrated Aesthetic Treatment System



Specifications subject to change without notice.
GentleMax, GentleFacial and Dynamic Cooling Device (DCD™) are trademarks and GentleYAG, GentleLASE and Candela are registered trademarks of Candela Corporation.
© 2007 Candela Corporation. All rights reserved. Printed in the USA. 07/07
0920-11-0054 Revision 06

Count On Candela To Get A Multi-Wavelength System Right

As the world leader in light-based aesthetic technologies, Candela continues to set the standard for innovation and performance. So it's no surprise that GentleMAX, our new integrated aesthetic system, offers more power and versatility than any other multi-wavelength device. With GentleMAX, you get a single consolidated system that delivers an incredibly diverse range of treatments—while also offering the option to add additional light-based technologies.

It Starts With Two Best-In-Class Lasers

GentleMAX brings together two gold standard lasers—GentleYAG and GentleLASE—to create one unbeatable system. GentleYAG, the fastest and most powerful Nd:YAG 1064 nm laser on the market, effectively treats all skin types, including tanned skin, for permanent hair reduction, leg and facial veins, wrinkles and skin tightening. GentleLASE, Candela's 755 nm laser, is the premier treatment for unwanted hair, pigmented and vascular lesions, and wrinkles. Individually, these lasers offer unrivaled flexibility and efficacy. Together, they deliver even more—allowing you to expand your treatment capabilities and enhance your practice with one stand-alone system.




Advanced Technology— True Power—Real Results

There are other multiple-wavelength systems on the market, but none come close to offering the efficacy, speed, and versatility of Candela's GentleMAX. When compared to the competition, it's easy to see the difference with Candela.

It Couldn't Be Any Easier

With Candela's exclusive touchscreen interface, choosing the ideal settings is as easy as 1-2-3. The GentleMax Guided Mode walks the user through the set-up process with a series of questions to help determine the optimal parameter settings based on individual patient criteria.

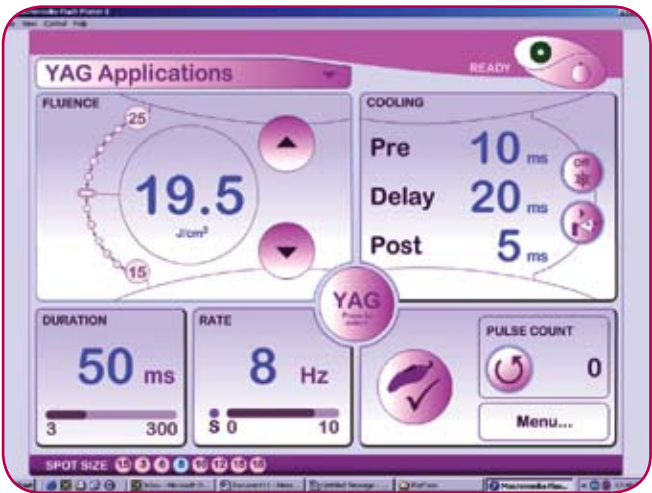

Your Choice of Epidermal Protection

With GentleMax you have your choice of skin cooling options—from chilled air technology to Candela's integrated and patented Dynamic Cooling Device that utilizes bursts of cryogen before and after the laser pulse to offer unparalleled patient comfort.



		and with 
hair reduction	•	•
permanent hair reduction	•	•
permanent hair reduction on all skin types	•	•
pseudofolliculitis barbae	•	•
wrinkle reduction	•	•
skin rejuvenation*	•	•
facial skin tightening	•	•
facial telangiectasia	•	•
PWS/hemangioma	•	•
cherry angiomas	•	•
venous lakes	•	•
leg veins	•	•
sun spots (solar lentigos)	•	•
age spots	•	•
ephiledes (freckles)	•	•
nevus of Ota/Ito	•	•
melasma	•	•
café-au-lait macules	•	•
blue nevi	•	•
poikiloderma of civatte	•	•
sebhorreic keratosis	•	•
non-facial skin tightening	•	•
warts	•	•
striae		•
rosacea		•
diffuse redness		•
spider angiomas		•
prolonged erythema		•
inflammatory acne		•
hypertrophic scars		•
skin resurfacing		•
skin tags		•
cutaneous lesions		•

*Skin rejuvenation by reduction of diffuse redness and pigmentation

The Option To Integrate Additional Technologies Offers Unprecedented Versatility

In addition to its 755 and 1064 wavelengths, GentleMAX can be expanded to meet your changing needs. Additional light-based technologies can be seamlessly integrated, creating a single system capable of delivering an incredibly wide range

of treatment applications. Other systems offer multi-wavelengths, but GentleMax is in a class by itself—offering you the power and performance you need to expand your practice and increase your profitability.



Integrating the SmoothPeel low-power erbium laser enables you to expand your aesthetic offering to include precise and predictable laser peels.

