

## BUILT WITH SAFETY IN MIND

- » Built-in audible indicators sound when each applicator reaches the optimal temperature, providing patients with updates on their treatment progress.
- » Evolve is designed with a patient call button to allow clients to easily pause their treatment and have clinic staff make adjustments if required.
- » Evolve's treatment screen enables physicians to visually monitor the thermal effect of each applicator at any point during the procedure.
- » Evolve's color blind RF technology ensures it can safely and effectively be used on all skin types.

Power Output:	Up to 50 W
Output Frequency:	1 MHz
Temperature Cut-Off:	Real time temperature monitoring for all active applicators. Operator selectable up to 43°C
Tailored Treatment Times:	5 minute increments allow procedures to be completely customized to specific patient needs



## EVOLVE ENABLES YOUR PRACTICE TO ADDRESS THE INDIVIDUAL TREATMENT NEEDS OF EACH BODY TYPE



ECTOMORPH

MESOMORPH

ENDOMORPH



INMODE IS COMMITTED TO BEING A LEADER YOU CAN TRUST.

PROVIDING UNCOMPARABLE CLINICAL, TECHNICAL & MARKETING SUPPORT.



TOLL FREE {855} 411 2639 | WWW.INMODEMD.COM

VSBEL0120

introducing



EVOLVE

by INMODE



TITE | TRIM | TONE  
HANDS-FREE TECHNOLOGY

All-in-One Body Treatments

## MAXIMIZE YOUR REVENUE & BOTTOM LINE

- » Evolve's **hands-free** design and programable technology allow your practice to reallocate physician and clinician time to other patients and treatments while Evolve is in use.
- » Evolve allows you to offer your patients skin remodelling, adipose tissue treatments, and muscle toning all with one platform. With Evolve you can provide multiple treatments to generate multiple revenue streams.



# THE NEXT GENERATION OF NON-INVASIVE, HANDS-FREE TECHNOLOGY DESIGNED TO REMODEL SKIN, TARGET ADIPOSE TISSUE AND TONE MUSCLES

According to The American Society for Aesthetic Plastic Surgery, in 2018 total non-invasive procedures grew at an impressive rate of +9% with over 3.3 million procedures performed<sup>1</sup>. Patients continue to pursue non-invasive treatments that improve skin appearance without the associated risks or downtime of invasive procedures.

Evolve is an innovative hands-free solution that maximizes patient results and provides the highest quality of aesthetic treatments available. This all-in-one platform, delivers multiple proprietary technologies on one system. The ultimate in medical grade technology, that ensures your practice capitalizes on the high demand body sculpting trend.

<sup>1</sup>-2018 The American Society of Aesthetic Plastic Surgery, Cosmetic (Aesthetic) Surgery National Data Bank.



## EVOLVE TITE HANDS-FREE SKIN REMODELING TREATMENT

Evolve Tite harnesses the power of radiofrequency energy to remodel skin and improve its appearance. The 8 hands-free applicators deliver uniform and volumetric heating to the skin & sub-dermal layer, to provide customized pain-free patient solutions.

Evolve's proprietary algorithm provides unprecedented safety. RF energy is automatically adjusted when each applicator senses it has reached the optimal thermal response.

Each applicator is designed with built-in thermal sensors which monitor skin surface temperature in real time allowing the device to reach an optimal thermal profile to deliver exceptional results.



DR. D. KIBA




## EVOLVE TRIM HANDS-FREE BODY & CELLULITE TREATMENT

Evolve Trim offers 6 hands-free applicators which deliver RF energy plus a negative tissue vacuum.

The patented RF energy and vacuum design deliver a uniform thermal effect to the deepest layers of subcutaneous fat to provide the maximum treatment.

Trim's 6 hands-free applicators enable your practice to simultaneously treat multiple body areas in one session.

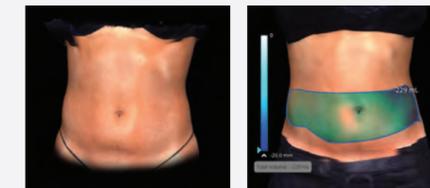

BEFORE & AFTER WITH QUANTIFICARE PHOTOGRAPHY  
DR. S. THEODOROU & DR. C. CHIA



## EVOLVE TONE HANDS-FREE MUSCLE TONING TREATMENT

Evolve Tone is a non-invasive muscle toning technology that works by emitting electromagnetic energy to stimulate involuntary muscle contractions.

Patients rest comfortably while Tone targets specific muscle groups to increase muscle strength and refine the look of muscles, without downtime.

Tone's 4 hands-free applicators allow for multiple body areas to be treated simultaneously, so your practice can deliver customized patient treatments. Provide your patients with the ultimate in tailored body treatments.



DR. H. RAMIREZ

